SAVE THESE INSTRUCTIONS

- Never operate the unit if it has a damaged plug, if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water. Refer to an Authorized Service Center for repair.
- DO NOT excessively bend or pull the cord. Unplugging the cord by grasping the plug – NOT THE CORD. Keep the cord away from heated surfaces.
- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Use this product only for its intended use as described in this manual. Only use attachments recommended by the manufacturer.
- Never use while drowsy.

- DO NOT submerge the unit in water or any other liquid.
- DO NOT place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
- DO NOT reach for a unit that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.
- The unit is designed for use in a 120 volt AC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT operate the unit in water or any other liquid.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING

- DO NOT use the unit while bathing.
- ALWAYS unplug this product immediately after using.
- DO NOT place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
- DO NOT reach for a unit that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.
- The unit is designed for use in a 120 volt AC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT submerge the unit in water or any other liquid.

DANGER
To reduce the risk of electrocution:

WARNING

To reduce the risk of burns, electrocution, fire, or injury:

- If you have had oral or periodontal surgery in the previous twelve months, contact your dentist prior to use.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if bleeding continues after two weeks of use.
- Be sure that only the water jet nozzle touches your teeth and gums. Damage to teeth and gums may occur if any other part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the Interplak® Dental Water Jet, and should always exercise particular care with any oral hygiene device.
- The Interplak® Dental Water Jet is a personal care unit and is intended for household use only.

MEDICAL WARNING


This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

OPERATING INSTRUCTIONS

Cleaning and Care

If, in the course of use, the valve (10) on the water tank begins to leak, this can be taken care of by cleaning the valve seat. To do this, push out the valve body using the nozzle tip, clean the valve seat and replace the valve.

- 1 Unit
- 2 Color-coded attachments:
A - Standard Tips
B - Subgingival Tips
C - Tongue Cleaner
- 3 Reservoir
- 4 Handle
- 5 Water stop button
- 6 Tip support
- 7 Rotary flow control
- 8 Drawer for tip storage
- 9 Power switch
- 10 Valve seat (inside tank)

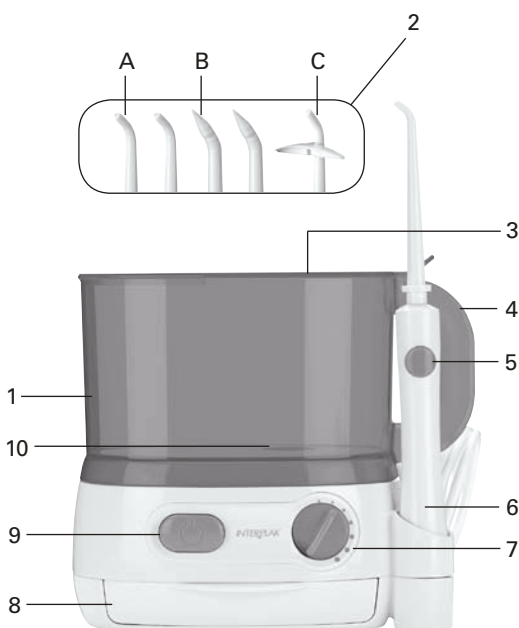Instructions For Use:

Remove the water reservoir (3). When removing reservoir, use one hand to securely hold base of the unit in place and, with other hand, lift tank straight up. Fill with lukewarm water. Reinstall the reservoir firmly onto the unit (1). Prior to initial use, the unit should be rinsed by letting a full reservoir of water pass through it. If the unit is to be used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for a long period of time, for reasons of hygiene it is advisable to follow the directions for initial use, letting a full reservoir of water pass through. To remove trapped air upon initial operation, turn the rotary flow control (7) all the way to the right and then to the left, three times.

Prior to filling the water reservoir, make sure water tip is firmly placed on the tip support and in the stand. The Water Jet is supplied with 5 color-coded attachments (2). Tips are located in convenient storage drawer (8).

Standard Tip (2pcs): removes trapped particles from between teeth.

Subgingival Tip (2pcs): uses a soft rubber tip to massage the gums and safely deliver solution into the gingival pockets below the gumline.

Tongue Cleaner (1pc): place the tongue cleaner at farthest end of tongue and move forward.

Remove the tip support (5) and install a tip onto the tip support. After unit has been switched on (9), the rotary flow control (7) can be used to regulate the water pressure. You should begin with the lowest setting so that your gums can become accustomed to the water jet.

Hold down the water-stop button (5) and, after switching on the unit, place the tip in your mouth and lean over the sink. When you release the water-stop button, the cleaning and massaging process will start. The direction of the water jet in your mouth can be controlled by turning the knurled ring on the slip-on nozzle. Keep your mouth slightly open and hold your elbow at mouth level, so that the water can run back into the sink. Direct the water jet toward your teeth, especially between your teeth, in order to rinse away stubborn remaining food.

This also massages your gums. Do not point the water jet directly toward your gums; point between your teeth to prevent damage to delicate gum areas.

Service Center:
Conair Corporation
Service Department
7475 N. Glen Harbor Blvd.
Glendale, AZ 85307
©2008 Conair Corporation, 150 Millford Road, East Windsor, NJ 08520
For more information regarding this and other Interplak® products, visit www.interplak.com

Limited 2-Year Warranty
Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service center listed below that is nearest you, together with your purchase receipt and \$3.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture. **ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.


Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Welcome to a New Level of At-Home Dental Care

The only water jet with built-in Microban® protection.

The Interplak® Water Jet enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the Water Jet can most effectively meet your specific oral needs.

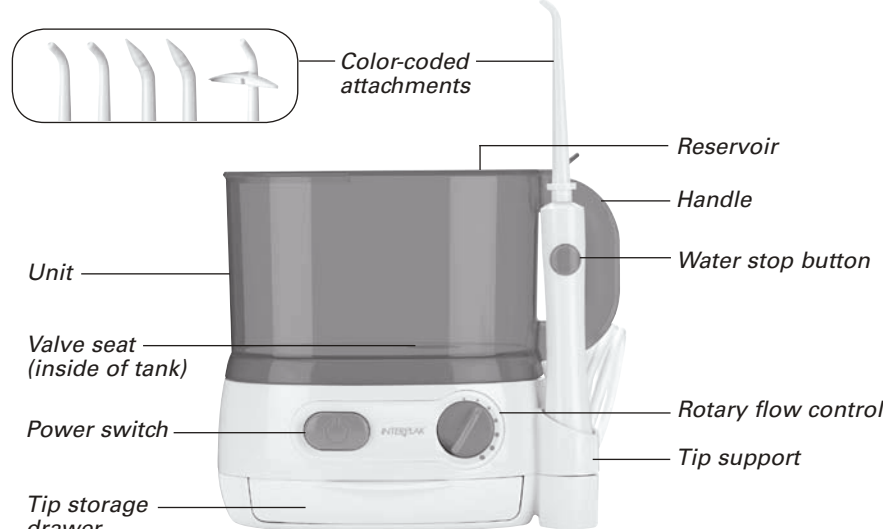DENTAL WATER JET

Instruction Manual • Model WJ10

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Questions or Comments? 1-800-334-4031
Visit us at our website: www.interplak.com

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